

Centripetal Forces & Hooke's Law




In today's lab, you will use the Welch Scientific Company's Centripetal Force Apparatus to verify the relationship between centripetal force, mass, velocity and radius. You may recall from class that for a mass moving in uniform circular motion, the centripetal force is given by:

$$F = m \frac{v^2}{R} \quad (\text{Equation 1})$$


where F is the centripetal force, m is the mass undergoing circular motion, R is the radius of the circular orbit and v is the tangential velocity of mass m .

A theoretical plot of the force as a function of mass versus $1/v^2$ would provide a straight line with a slope

which is equal to $(1/FR)$. For this portion of the experiment, both F and R would be constant: the force is equal to the amount of force required to hold the pendulum over the desired location and this occurs at the same radius each time. The plot of this would appear as shown.

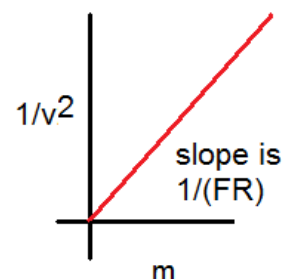

You will explore this for several different values of m in this lab, namely the mass of the pendulum. I will call this mass m_p . You will (by screwing masses onto the top of the pendulum) also investigate it for m_p , $m_p + 50g$, $m_p + 100g$, $m_p + 150g$. Experimental note: when you run your experiments, the upper most screw must be tight and the correct position for the pendulum is over top of the indicator bar with the spring, string, bar, and stand **making a square**. If these two things are not done, you will have larger experimental errors.

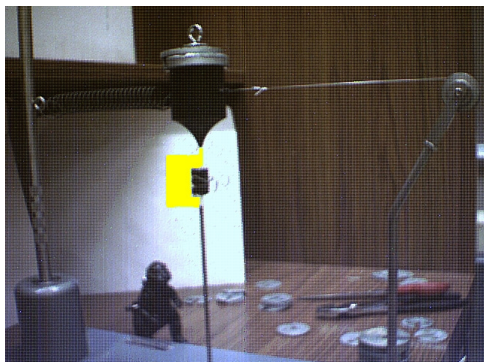

We could also vary the radius of orbit but this proves to be significantly more difficult so we will keep R fixed in this lab. The image to the left shows the apparatus.

There is a small tab of paper attached to the upward sticking rod. When the apparatus is spinning at the correct velocity so that the pendulum mass is aligned over the metal indicator, and the mass is spinning, you will hear a click as the mass strikes the paper tab as shown in the image to the left.

Procedure: For each of the 4 masses, you should find the time T_n that it takes the system to spin through n revolutions. The period (T) is then given by

$$T = \frac{T_n}{n}$$

Here, n should be at least 20, but you can do more if you are so inclined. T is the amount of time for one revolution. Next, you should measure the distance from the center to the tab of paper. This will be R . Be sure to convert this measurement into meters. Now, in order to find out how much force is required to move the mass to the tab of paper, use the weight hanger together with the special attachment which I have provided (you can see this in the second figure). Add weight as shown in the first image until the pendulum is aligned with the tab of paper.

A question now: what are you measuring when you measure this particular force which is given by:

$$F = mg?$$

The answer is that you are measuring how much force the spring is exerting on the pendulum while it is swinging trying to pull the mass backward towards the center. Of course, you're not going to want to be spinning the thing until you remove that paperclip attachment that I've made for you!

Now the next thing is this: How do you know what the tangential velocity is that you are measuring? And here is the answer: in a period T , the mass travels through a distance given by:

$$s = 2\pi R$$

Thus you are able to find the tangential velocity pretty easily:

$$v = \frac{2\pi R}{T}$$

Important experimental notes

Here are important notes to reduce your errors: when the spring is removed from the pendulum, the pendulum should hang right over the rod that is sticking up. You will need to adjust the screw on top to make sure this is the case. Failure to do this will produce larger errors. After reconnecting the spring to the pendulum, you will want to balance the system by moving the adjustable mass on top. The system is balanced when it rotates without jumping around. You will also want to make sure the system is leveled by adjusting the legs on the board.

Analysis: The centripetal force is related to the mass, velocity and radius by the first equation above. For each of your measurements, calculate the experimental centripetal force from Equation 1. Compare these calculations to your measured centripetal force by using percent error. You may consider that the correct value for the centripetal force is obtained by the masses attached to the hanger.

In this lab, you will also want to plot a graph of m vs. $(1/v^2)$. Looking at the theoretical equation for the centripetal force, you would expect the slope of this graph to be equal to

$$\text{slope} = \frac{1}{RF}$$

This is how this comes about: you measure m , v and R . Using the first equation, you can then obtain:

$$F = m \frac{v^2}{R} \Rightarrow \frac{1}{v^2} = \frac{m}{RF}$$

In this experiment, you are varying m and observing the effect that this has upon v . Now let's define a new variable, y such that:

$$y \equiv \frac{1}{v^2}$$

In terms of this variable, the equation looks like:

$$y = \left[\frac{1}{RF} \right] m = [\text{slope}] m$$

So if you plot a graph of m vs. y , the slope of this graph will be

$$\text{slope} = \frac{1}{RF}$$

You might ask: why worry about making this linear? One easy answer is this makes the fitting process simpler. Now once the slope has been obtained, you can determine the force from:

$$F = \frac{1}{R[\text{slope}]}$$

Compare your slope to the theoretical value of F that is measured by hanging weights over the end of it as shown in the second figure. In principle, you should obtain fairly good agreement here. Discuss in your write-up where experimental errors may have entered into the results.

Part II: Hooke's Law

In this part of the lab, do not stretch the silver springs beyond 10 cm. You should wear eye protection during this part of the lab also.

As you have seen, springs are very useful for measurement of force. You may recall from class that to first order, the force required to stretch a spring is related to the distance the spring is stretched by Hooke's law:

$$|\vec{F}_{\text{spring}}| = k |\Delta \vec{x}| \text{ (Equation 2)}$$

where Δx is the displacement from equilibrium and k is the spring constant. I would like for you to learn how to quickly measure the spring constant of springs. Although you will not be using this in today's lab, it is an important technique that you will use in at least one later lab. You will need to construct several setups for the measurement of the spring constant for 3 different springs (1 large spring, 1 small spring and two small springs in parallel). This is so you can compare one of your other springs to this system and see if the spring constants for springs in parallel add or if the inverses add. **You must answer the following question based upon your data:** do springs in parallel add as: $\frac{1}{k_{\text{eff}}} = \frac{1}{k_1} + \frac{1}{k_2}$ or do they add as: $k_{\text{eff}} = k_1 + k_2$?

In today's lab, I want you to measure the spring constant of the large copper spring, the small silver spring and the slinky. I want you to confirm how spring constants add

only for the silver spring. When you measure the spring constant of the slinky, do it in the following way: find the mass of the slinky, count the number of turns on the slinky and then calculate the mass per turn by dividing the two. With a paper board, you can then hang 5 turns, 10 turns , (to get the mass hanging, you multiply the mass per turn by the number of turns hanging) measuring the hanging length each time. Plot your results on a graph with force on the x-axis and length on the y axis and provide a qualitative discussion of the results.

You should measure the spring constant by use of a mass hanger and application of mass. You will need to make a data table containing m , F and x in SI units for the 3 springs. For 10 different masses, record the displacement of the spring and then for each spring, plot a graph of x vs. F . Other than performing a good experiment here and determining your spring constants, no further analysis is required for part II except for the determination that I mentioned above. However, you will most certainly want to write about this in your lab writeup. You will also want to remember how to do this for future reference.